






Improving Quality of Life



X-Ray Grid

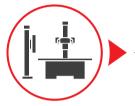

Radiography System



Software



Computed Radiography



Processor



X-Ray Accessories



RF Treatment





Economic Digital Radiography System



Clear Vision





Operation Software

ExamVue software presents a comprehensive set of X-ray image analysis tools through the use of an fast and easy-to-use interface. ExamVue offers high performance image processing and measurement functions.

X-ray Generator

High frequency generator ensures minimum patient dose, excellent reproducibility, and superior image contrast.

Flat panel detector

Flat panel detector for digital radiographic applications delivers high sensitivity, high resolution, low noise and superb image quality. Portable flat panel detector provides unlimited clinical applications.

Tube stand

Ergonomically designed tube stand provides the operator with faster and easier movement for table and stand examinations.

Bucky stand

Wide range of vertical travel ($400 \sim 1,800$ mm) allows exposures to be carried out from the cervical spine to the lower extremities.

Bucky table

The robust table design allows the examination of obese patients, with a table weight capacity of up to 150kg. Due to the footswitch integrated into the table base the operator can care the patient safely without any unexpected movement of the tabletop.







The system can be used as a dual 17" X 17" detector system or a single 14" X 17" detector system.

Specifications of the system

Generator	DRE140	DRE140 DRE150	
Max. Output	40kW	50kW	
mA range	10 ~ 500mA	10 ~ 630mA	
kV range	40 ~ 125kV	40 ~ 150kV	
mAs range	0.4 ~ 630mAs	0.1 ~ 630mAs	
Power supply	1Ø 220VAC 50/60Hz	3Ø 380VAC 50/60Hz	
Exposure time	0.004 ~ 10sec	0.001 ~ 6.3sec	
No. of APR	216	1,024	
Option	AEC	AEC/ High speed rotor/ DAP	

X-ray tube	DRE140 DRE150			
Focal spot	1.0/2.0mm	0.6/1.2mm		
Anode heat capacity	100kJ (140kHU)	210kJ (300kHU)		
Target angle	16°	12°		
Voltage range	40 - 125kV	40 - 150kV		
Weight	16kg	18kg		
Permanent filtration	0.9mmAl	0.9mmAl		

Specifications of the system

Flat panel detector	DRE140/150	40/150 DRE140/150		
Scintillator	CsI	CsI		
Pixel pitch	143µm 143µm			
Resolution	3.7lp/mm	3.7lp/mm		
MTF	36% (2.0lp/mm, 70kVp, 1x1)	36% (2.0lp/mm, 7 0kVp, 1x1)		
DQE	70%	70%		
Active Area	43 x 43.9cm	35 x 43cm		
Pixel matrix	3,008 x 3,072	2,448 x 2,984		
AD conversion	14 bits	16 bits		
Energy range	40 ~ 150kVp	40 ~ 150kVp		
Power supply	1Ø 100 ~ 240VAC 50/60Hz	1Ø 100 ~ 240VAC 50/60Hz		
Dimension	512 x 495 x 43mm	460 x 384 x 15mm		
Weight	16kg	3kg		

Collimator	DRE140/150
Power supply	1Ø 24VAC 50/60Hz
Max. kVp shield	150kV
Max. field size	35 x 35cm at SID 65cm
Inherent filtration	1.2mmAl
Lamp luminosity	Over 160Lux
Lamp time	30sec
Field size control	Manual

Specifications of the system

Tube stand

Vertical travel distance	1,400mm	
Transverse travel distance	260mm	
Min. Tube height (from floor)	550mm	
Max. Tube height (from floor)	1,950mm	
Column rotation	±90°	~~
Tube rotation	±180°	
Longitudinal travel distance	2,400 ~ 3,000mm (floor ceiling fixed type) 1,500mm (floor fixed type)	
Power supply	24VDC, 3A	
Weight	200kg (with rails)	
Application	DRE140/150	

Bucky stand

Vertical travel distance	1,400mm	
Min. Bucky height (from floor)	400mm	
Max. Bucky height (from floor)	1,800mm	
Power supply	1Ø 16 ~ 24VAC 50/60Hz	
Weight	120kg	
Application	DRE140/150	

Bucky table

Table top size	2,000 x 745mm
Table top height	670mm
Table top material	Acryl
Longitudinal Travel distance	±280mm
Transverse travel distance	±140mm
Bucky travel distance	±230mm
Max. patient weight	Max. 150kg
Power supply	1Ø 220VAC 50/60Hz
Weight	160kg
Application	DRE140/150

Types of system



Clear Vision DRE 140/150

Single Detector (14" X 17" or 17" X 17")

Bucky Stand

+

Tube Stand system

Configuration

Bucky stand, Tube stand(floor mounted or floor ceiling),

Generator, PC&Monitor, Acquire software, Flat Panel Detector, Collimator

Clear Vision DRE 140/150

Single Detector (14" X 17" or 17" X 17")

Tube Stand

+

Bucky Table(4way) system

Configuration

4Way Bucky table, Tube stand(floor mounted or floor ceiling),

Generator, PC&Monitor, Acquire software, Flat Panel Detector, Collimator

Clear Vision DRE 140/150

Single Detector (14" X 17" or 17" X 17")

Dual Detector (14" X 17" or 17" X 17")

Bucky Stand

+

Tube Stand

+

Bucky Table(4way) system

Configuration

4Way Bucky table, Bucky stand, Tube stand(floor mounted or floor ceiling), Generator, PC&Monitor, Acquire software, Flat Panel Detector, Collimator

Clear Vision DRE series


Economic Digital Radiography System

Configuration of the system

Components	DRE140 with fixed dual FPD	DRE140 with portable FPD (★ with CR/film)	DRE150 with fixed dual FPD	DRE150 with portable FPD (★ with CR/film)
40kW Generator	~	~ ~	_	_
50kW Generator	-	-	~	~ ~
17" x 17" Flat panel detector	~	-	~	-
14" x 17" Flat panel detector	-	~	-	~
150kV X-ray tube	-	-	~	~ ~
125kV X-ray tube	~	~	_	- -
Basic Configuration	★ Tube stand, Bucky stand, Bucky table, PC & monitor, Acquisition Software			

- ★ upgradable to floor mounted Tube stand
- ★ for international sales



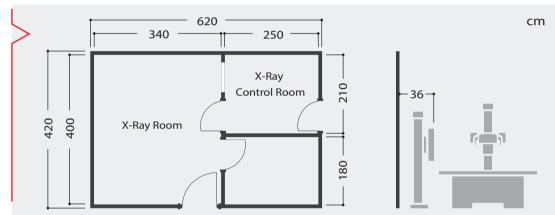

ExamVue Software

for Medical







ExamVue for Med

Imaging and Annotation Tools

(Window level, Zoom, Pan, Crop, Invert, Rotation, Angle, Label, Text)

Powerful User Features

(Image stitching, Image import / export, Image comparison, Multi view, Image expansion, Full screen, Auto collimation, Reject/Accept)

Smart Window Level

DICOM3.0 Compatibility

Post Image processing mode for each body part

Changeable Skin color (Green Tea, Strawberry, Mint, Chocolate, Banana, Vanilla)

VK FORCE SOLUTIONS SDN. BHD. (1288081-H) No. 933, Block A3, Jalan PJS 8/9, Leisure Commerce Square, 46150 Petaling Jaya, Selangor Malaysia.

Tel: +60378654913

Sales Hotline: +6018-295 8292 | +6018-235 8292 www.vkforcesolutions.com

Email: sales@vkforcesolutions.com

Registered:









